

# RÜSCH ENDOTRACHEAL TUBES

The optimum tube for all applications. At RÜSCH, this demand is met thanks to an extraordinarily extensive and versatile range of tubes. All RÜSCH tubes feature an outstanding I.D. to O.D. ratio. Our high quality standards are reflected in the high-grade materials we use, which are tested according to the highest medical standards.

To meet the special requirements of paediatric patients, Teleflex's RÜSCH brand supplies a large variety of tracheal tubes in smaller sizes. Tracheal tubes for paediatric care are marked with black tip to ensure safe tracheal positioning.

CHILD AGE	TUBE I.D.	BLACK TIP LENGTH
1.	2.0–3.5 mm	20 mm
2.–5.	4.0–5.0 mm	30/40 mm
6.–14.	5.5–6.5 mm	40 mm

## TRACHEAL TUBES:

1. Before using tracheal tubes with the connector pushed on halfway, it must be ensured that the connector is pushed into the tube shaft as far as it will go.
2. To improve the tight hold, it is advisable to clean with ethanol, both the contact surface of a connector that has been loosened once or removed and the corresponding tube shaft, before reinserting it.
3. This particularly applies to those cases in which the connector was completely removed first (e.g. for fiberoptic insertions), or was pushed on again after the tube had been shortened.

## SUPER SAFETYCLEAR



### CLEAR TRACHEAL TUBE MADE OF PVC, WITH LOW-PRESSURE CUFF

nasal/oral

sizes: I.D. 2.5–10.0 mm

semi-seated connector, valve for Luer and Luer-lock syringes, cupped atraumatic tip, continuous X-ray marker, black position indicator for correct tube placement, blue pilot balloon, graduated

- latex-free
- sterile



### SUPER SAFETYCLEAR

RÜSCH

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY	
112480, Magill	2.5 mm	4.0 mm	8.0 mm	160 mm	10	
	3.0 mm	5.0 mm	8.0 mm	170 mm		
	3.5 mm	5.3 mm	8.0 mm	190 mm		
	4.0 mm	6.0 mm	10.5 mm	220 mm		
	4.5 mm	6.3 mm	10.5 mm	230 mm		
	112482, Murphy, 1 eye	5.0 mm	6.7 mm	13.0 mm		250 mm
		5.5 mm	7.3 mm	16.5 mm		280 mm
		6.0 mm	8.0 mm	18.5 mm		290 mm
		6.5 mm	8.7 mm	20.5 mm		300 mm
		7.0 mm	9.3 mm	24.0 mm		320 mm
7.5 mm		10.0 mm	26.0 mm	330 mm		
8.0 mm		10.7 mm	26.0 mm	340 mm		
8.5 mm		11.3 mm	28.0 mm	345 mm		
9.0 mm		12.0 mm	28.0 mm	350 mm		
9.5 mm		12.7 mm	29.0 mm	350 mm		
10.0 mm	13.3 mm	29.0 mm	350 mm			

\* Length without connector

## FLEXISET



### SUPER SAFETYCLEAR TRACHEAL TUBE WITH MURPHY EYE AND FLEXISLIP INTUBATION STYLET

sizes: I.D. 5.0–9.0 mm

To order each product separately

- Super SafetyClear tracheal tube only, Ref. 112482
- FlexiSlip intubation stylet only, Ref. 503700  
(please refer to page 56)

#### FLEXISET

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
112481	5.0 mm	6.7 mm	13.0 mm	250 mm	10
	5.5 mm	7.3 mm	16.5 mm	280 mm	
	6.0 mm	8.0 mm	18.5 mm	290 mm	
	6.5 mm	8.7 mm	20.5 mm	300 mm	
	7.0 mm	9.3 mm	24.0 mm	320 mm	
	7.5 mm	10.0 mm	26.0 mm	330 mm	
	8.0 mm	10.7 mm	26.0 mm	340 mm	
	8.5 mm	11.3 mm	28.0 mm	345 mm	
9.0 mm	12.0 mm	28.0 mm	350 mm		

\* Length without connector



## SAFETYCLEAR PLUS

### CLEAR TRACHEAL TUBE MADE OF PVC WITH HIGH-VOLUME LOW-PRESSURE SUPER PLUS CUFF nasal/oral

sizes: I.D. 5.0–10.0 mm

semi-seated connector, valve for Luer and Luer-lock syringes, cupped atraumatic tip, continuous X-ray marker, black position indicator for correct tube placement, blue pilot balloon, graduated

- latex-free
- sterile



#### SAFETYCLEAR PLUS

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
112080, Magill	5.0 mm	6.7 mm	20 mm	250 mm	10
	5.5 mm	7.3 mm	20 mm	280 mm	
	6.0 mm	8.0 mm	24 mm	290 mm	
	6.5 mm	8.7 mm	24 mm	300 mm	
	7.0 mm	9.3 mm	28 mm	320 mm	
	7.5 mm	10.0 mm	28 mm	330 mm	
	8.0 mm	10.7 mm	30 mm	340 mm	
	8.5 mm	11.3 mm	30 mm	345 mm	
112082, Murphy, 1 eye	9.0 mm	12.0 mm	34 mm	350 mm	
	9.5 mm	12.7 mm	34 mm	350 mm	
	10.0 mm	13.3 mm	34 mm	350 mm	

\* Length without connector

# SAFETYCLEAR




## CLEAR TRACHEAL TUBE MADE OF PVC WITHOUT CUFF

nasal/oral  
sizes: I.D. 2.0–10.0 mm

semi-seated connector, cupped atraumatic tip, continuous X-ray marker, graduated

- latex-free
- sterile

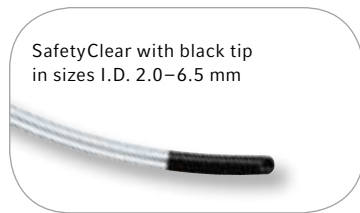

SafetyClear with black tip  
in sizes I.D. 2.0–6.5 mm



### SAFETYCLEAR

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	LENGTH BLACK TIP	QTY
100380, Magill sizes: I.D. 2.0–10.0 mm	2.0 mm	3.0 mm	150 mm	20 mm	10
	2.5 mm	3.3 mm	160 mm	20 mm	
	3.0 mm	4.0 mm	170 mm	20 mm	
	3.5 mm	4.7 mm	190 mm	20 mm	
	4.0 mm	5.3 mm	210 mm	30 mm	
	4.5 mm	6.0 mm	230 mm	30 mm	
100382, Murphy, 1 eye sizes: I.D. 2.0–7.0 mm	5.0 mm	6.7 mm	250 mm	40 mm	
	5.5 mm	7.3 mm	280 mm	40 mm	
	6.0 mm	8.0 mm	290 mm	40 mm	
	6.5 mm	8.7 mm	300 mm	40 mm	
	7.0 mm	9.3 mm	320 mm	–	
	7.5 mm	10.0 mm	330 mm	–	
	8.0 mm	10.7 mm	340 mm	–	
	8.5 mm	11.3 mm	345 mm	–	
	9.0 mm	12.0 mm	350 mm	–	
	9.5 mm	12.7 mm	350 mm	–	
	10.0 mm	13.3 mm	350 mm	–	

\* Length without connector

## ISIS ETT

### ENDOTRACHEAL TUBE WITH INTEGRATED PORT FOR OPTIONAL SUBGLOTTIC SECRETION SUCTIONING

The ideal tube to provide the best of care to all your patients, up to 30 days.

nasal/oral 1 eye (Murphy)

Sizes: I.D. 6.0–9.0 mm

This convertible endotracheal tube with separate suction line allows for subglottic secretion removal on demand. This frees the clinician from the burden of deciding for an expensive tube at the time of intubation or to extubate/re-intubate a patient to provide a VAP prevention method of subglottic secretion suctioning. Thanks to the ISIS flexible design, subglottic suctioning can be provided to any patient whenever necessary.

### SUCTION LINE FOR ISIS ENDOTRACHEAL TUBE

When connected to the ISIS endotracheal tube, the optional suction line provides access from the top of the ETT cuff for easy & safe suctioning of the subglottic secretions that spill over the cuff. Available separately.

1 size fits all ISIS tubes



#### TELEFLEX ISIS ETT

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
112662	6.0 mm	9.2 mm	20.5 mm	290 mm	10
	6.5 mm	9.9 mm	24 mm	300 mm	
	7.0 mm	10.6 mm	26 mm	320 mm	
	7.5 mm	11.3 mm	26 mm	330 mm	
	8.0 mm	12.0 mm	28 mm	340 mm	
	8.5 mm	12.6 mm	28 mm	345 mm	
	9.0 mm	13.2 mm	29 mm	350 mm	

\* Length without connector

Endotracheal tube with integrated port for optional subglottic secretion suctioning.



#### SUCTION LINE

REF.	QTY
502700	5

# SUPER SAFETY SILK



**TRACHEAL TUBE  
MADE OF DEHP-FREE PVC  
WITH LOW PRESSURE  
HIGH VOLUME CUFF**

nasal/oral  
sizes: I.D. 3.0–10.0 mm

tissue-friendly thermosensitive material, especially advantageous for long-term intubation, soft atraumatic tip, continuous X-ray marker, semi-seated connector, low pressure sealing, Safety Silk surface, black ring for intubation depth control

- latex-free
- sterile



## SUPER SAFETY SILK

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH	QTY
112682, Murphy, 1 eye sizes: I.D. 3.0–10.0 mm	3.0 mm	5.3 mm	8.0 mm	170 mm	10
	3.5 mm	6.0 mm	10.5 mm	190 mm	
	4.0 mm	6.3 mm	10.5 mm	220 mm	
	4.5 mm	6.7 mm	13.0 mm	230 mm	
112680, Magill sizes: I.D. 5.0–10.0 mm	5.0 mm	7.3 mm	16.5 mm	250 mm	
	5.5 mm	8.0 mm	18.5 mm	280 mm	
	6.0 mm	8.7 mm	20.5 mm	290 mm	
	6.5 mm	9.3 mm	24.0 mm	300 mm	
	7.0 mm	10.0 mm	26.0 mm	320 mm	
	7.5 mm	10.7 mm	26.0 mm	330 mm	
	8.0 mm	11.6 mm	28.0 mm	340 mm	
8.5 mm	12.4 mm	28.0 mm	345 mm		
9.0 mm	13.0 mm	28.0 mm	350 mm		
9.5 mm	14.2 mm	29.0 mm	350 mm		
10.0 mm	15.0 mm	29.0 mm	350 mm		

## PEDI SAFETY SILK



### CLEAR TRACHEAL TUBE MADE OF SOFT PVC

part of the Safety Silk family,  
nasal/oral (Magill)

sizes: I.D. 2.0–5.5 mm

semi-seated connector, cupped atrau-  
matic black tip, continuous X-ray  
marker, soft version, graduated

- latex-free
- sterile




### PEDI SAFETY SILK

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	BLACK TIP LENGTH	QTY
100480	2.0 mm	3.3 mm	140 mm	20 mm	10
	2.5 mm	4.0 mm	160 mm	20 mm	
	3.0 mm	4.7 mm	180 mm	20 mm	
	3.5 mm	5.3 mm	200 mm	20 mm	
	4.0 mm	6.0 mm	220 mm	30 mm	
	4.5 mm	6.7 mm	240 mm	30 mm	
	5.0 mm	7.3 mm	260 mm	40 mm	
	5.5 mm	8.0 mm	280 mm	40 mm	

\* Length without connector

## NASAL SAFETY SILK

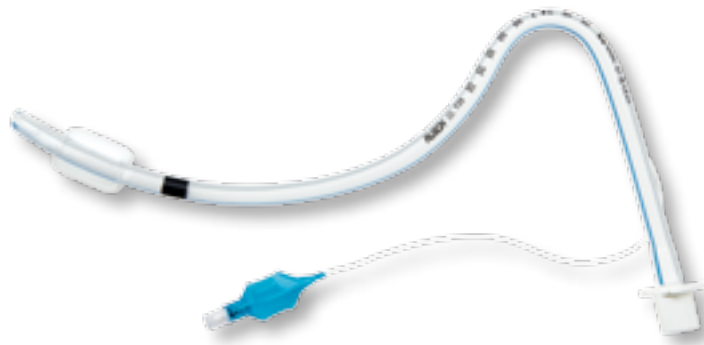

### CLEAR SOFT ANATOMICALLY SHAPED TUBE

MADE OF DEHP-FREE PVC

nasal only, 1 eye (Murphy)

sizes: I.D. 3.0–8.0 mm

extremely soft & tissue-friendly PVC,  
bonded connector, kink-resistant,  
cupped atraumatic tip, continuous blue  
X-ray marker, graduated, with low-  
pressure cuff, blue pilot balloon, valve  
for Luer and Luer-lock syringes

- latex-free
- sterile




### NASAL SAFETY SILK

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH	QTY
111782	3.0 mm	5.3 mm	11.0 mm	240 mm	10
	3.5 mm	6.0 mm	11.0 mm	271 mm	
	4.0 mm	6.3 mm	13.0 mm	303 mm	
	4.5 mm	6.7 mm	13.0 mm	338 mm	
	5.0 mm	7.3 mm	16.0 mm	358 mm	
	5.5 mm	8.0 mm	18.0 mm	368 mm	
	6.0 mm	8.7 mm	20.0 mm	415 mm	
	6.5 mm	9.3 mm	24.0 mm	420 mm	
	7.0 mm	10.0 mm	26.0 mm	430 mm	
	7.5 mm	10.7 mm	26.0 mm	445 mm	
8.0 mm	11.3 mm	28.0 mm	455 mm		

# AGT TRACHEAL TUBE ORAL

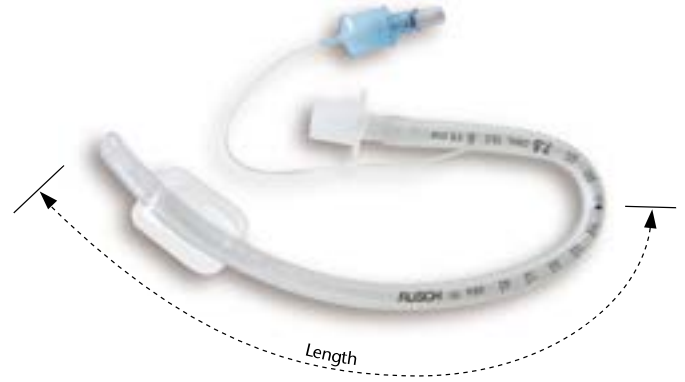

## CLEAR ANATOMICALLY SHAPED AGT TUBE MADE OF PVC WITH LOW-PRESSURE CUFF


oral, 1 eye  
 sizes: I.D. 3.5–9.0 mm

semi-seated connector, cupped atraumatic tip, continuous X-ray marker, graduated, low-pressure cuff, blue pilot balloon, valve for Luer and Luer-lock syringes

- latex-free
- sterile



AGT TRACHEAL TUBE, ORAL					RÜSCH
REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
111780	3.5 mm	5.3 mm	11 mm	120 mm	10
	4.0 mm	6.0 mm	11 mm	130 mm	
	4.5 mm	6.3 mm	13 mm	145 mm	
	5.0 mm	6.7 mm	13 mm	155 mm	
	5.5 mm	7.3 mm	16 mm	165 mm	
	6.0 mm	8.0 mm	18 mm	180 mm	
	6.5 mm	8.7 mm	20 mm	190 mm	
	7.0 mm	9.3 mm	24 mm	200 mm	
	7.5 mm	10.0 mm	26 mm	210 mm	
	8.0 mm	10.7 mm	26 mm	220 mm	
	8.5 mm	11.3 mm	28 mm	230 mm	
9.0 mm	12.0 mm	28 mm	240 mm		

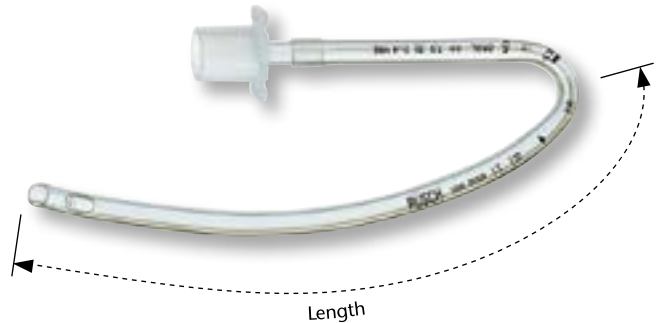

## CLEAR ANATOMICALLY SHAPED AGT TUBE MADE OF PVC WITHOUT CUFF

oral, 2 eyes  
 sizes: I.D. 3.0–7.0 mm

semi-seated connector, cupped atraumatic tip, continuous X-ray marker, graduated

- latex-free
- sterile

AGT TRACHEAL TUBE, ORAL					RÜSCH	
REF.	ORDER SIZE/I.D.	O.D.	O.D. (CH.)	LENGTH*	RECOMMENDED SIZE FOR	QTY
100180	3.0 mm	4.0 mm	Ch. 12	105 mm	newborn babies	10
	3.5 mm	4.7 mm	Ch. 14	120 mm	1–6 months	
	4.0 mm	5.3 mm	Ch. 16	130 mm	0.5–1 year	
	4.5 mm	6.0 mm	Ch. 18	145 mm	1–2 years	
	5.0 mm	6.7 mm	Ch. 20	155 mm	2–4 years	
	5.5 mm	7.3 mm	Ch. 22	165 mm	4–6 years	
	6.0 mm	8.0 mm	Ch. 24	180 mm	6–8 years	
	6.5 mm	8.7 mm	Ch. 26	190 mm	8–9 years	
7.0 mm	9.3 mm	Ch. 28	200 mm	over 9		

\* For length see picture


## AGT TRACHEAL TUBE NASAL

### CLEAR ANATOMICALLY SHAPED AGT TUBE MADE OF PVC WITH LOW-PRESSURE CUFF

nasal only, 1 eye  
sizes: I.D. 3.5–8.0 mm

semi-seated connector, cupped atraumatic tip, continuous X-ray marker, graduated, low-pressure cuff, blue pilot balloon, valve for Luer and Luer-lock syringes


- latex-free
- sterile

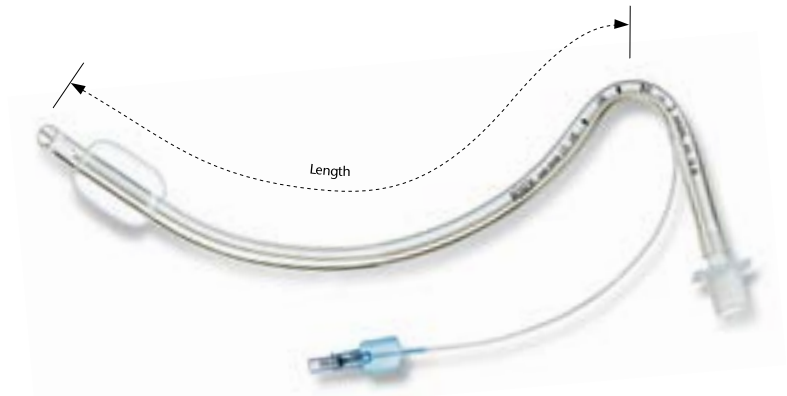

### AGT TRACHEAL TUBE, NASAL

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
111781	3.5 mm	5.3 mm	11 mm	160 mm	10
	4.0 mm	6.0 mm	11 mm	185 mm	
	4.5 mm	6.3 mm	13 mm	220 mm	
	5.0 mm	6.7 mm	13 mm	235 mm	
	5.5 mm	7.3 mm	16 mm	245 mm	
	6.0 mm	8.0 mm	18 mm	260 mm	
	6.5 mm	8.7 mm	20 mm	270 mm	
	7.0 mm	9.3 mm	24 mm	280 mm	
	7.5 mm	10.0 mm	26 mm	290 mm	
8.0 mm	10.7 mm	26 mm	300 mm		

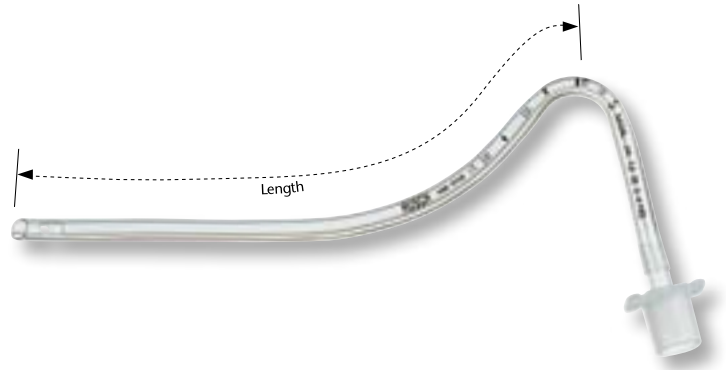

### CLEAR ANATOMICALLY SHAPED AGT TUBE MADE OF PVC WITHOUT CUFF

nasal only, 2 eyes  
sizes: I.D. 3.0–7.0 mm

semi-seated connector, cupped atraumatic tip, continuous X-ray marker, graduated

- latex-free
- sterile

### AGT TRACHEAL TUBE, NASAL

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	O.D. (CH.)	LENGTH*	RECOMMENDED SIZE FOR	QTY
100181	3.0 mm	4.0 mm	Ch. 12	115 mm	newborn babies	10
	3.5 mm	4.7 mm	Ch. 14	140 mm	1–6 months	
	4.0 mm	5.3 mm	Ch. 16	155 mm	0.5–1 year	
	4.5 mm	6.0 mm	Ch. 18	180 mm	1–2 years	
	5.0 mm	6.7 mm	Ch. 20	190 mm	2–4 years	
	5.5 mm	7.3 mm	Ch. 22	195 mm	4–6 years	
	6.0 mm	8.0 mm	Ch. 24	205 mm	6–8 years	
	6.5 mm	8.7 mm	Ch. 26	210 mm	8–9 years	
7.0 mm	9.3 mm	Ch. 28	215 mm	over 9		

\* For length see picture



# SILKOCLEAR



## CLEAR TRACHEAL TUBE MADE OF SILICONE WITH PREFORMED CUFF

nasal/oral, 1 eye (Murphy)  
sizes: I.D. 5.0–9.0 mm

fixed connector, valve for Luer and Luer-lock adapter syringes, cupped atraumatic tip, continuous X-ray marker, black position indicator for correct tube placement, transparent pilot balloon, fluted inner lumen, graduated

- latex-free
- sterile

### SILKOCLEAR

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
105102	5.0 mm	7.3 mm	8.5 mm	240 mm	5
	5.5 mm	8.0 mm	12.0 mm	270 mm	
	6.0 mm	8.7 mm	12.0 mm	280 mm	
	6.5 mm	9.3 mm	13.5 mm	290 mm	
	7.0 mm	10.0 mm	13.5 mm	300 mm	
	7.5 mm	10.7 mm	13.5 mm	310 mm	
	8.0 mm	11.3 mm	16.0 mm	320 mm	
	8.5 mm	12.0 mm	16.0 mm	330 mm	
	9.0 mm	12.7 mm	16.0 mm	340 mm	



## CLEAR TRACHEAL TUBE MADE OF SILICONE WITHOUT CUFF

nasal/oral (Magill)  
sizes: I.D. 2.0–4.5 mm

semi-seated connector, black tip,  
graduated

- latex-free
- sterile

### SILKOCLEAR

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	BLACK TIP LENGTH	QTY
100280	2.0 mm	3.3 mm	120 mm	20 mm	10
	2.5 mm	4.0 mm	140 mm	20 mm	
	3.0 mm	4.7 mm	160 mm	20 mm	
	3.5 mm	5.3 mm	180 mm	20 mm	
	4.0 mm	6.0 mm	200 mm	30 mm	
	4.5 mm	6.7 mm	220 mm	30 mm	

\* Length without connector

**Please note:**

SILKOSPRAY is not suitable for devices made of silicone. If you require a lubricant, we recommend our TrachJell (page 59) or commercially available gels.

# RÜSCHFLEX

## ARMOURED PREFORMED TRACHEAL TUBE, CURVED, MADE OF PVC WITH LOW-PRESSURE CUFF

nasal/oral

sizes: I.D. 3.5–10.0 mm

fixed white connector, blue pilot balloon, valve for Luer- and Luer-lock syringes, radiopaque, graduated, also available with integrated FlexiSlip stylet (can be ordered separately, Ref. 503700. see page 56 for more details)

- latex-free
- sterile



*Magill, Ref. 104201  
Size: I.D. 3.5–10.0 mm*



*Magill, with FlexiSlip,  
Ref. 104203  
Size: I.D. 5.0–9.0 mm*



## RÜSCHFLEX

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH* MAGILL 104201/104203	LENGTH* MURPHY 104202/104204	QTY
104201, Magill sizes: I.D. 3.5–10.0 mm	3.5 mm	5.8 mm	228 mm	229 mm	5
	4.0 mm	6.3 mm	235 mm	237 mm	
104203, Magill with FlexiSlip sizes: I.D. 5.0–9.0 mm	4.5 mm	6.8 mm	256 mm	259 mm	
	5.0 mm	7.3 mm	319 mm	321 mm	
	5.5 mm	7.8 mm	319 mm	322 mm	
	6.0 mm	8.4 mm	326 mm	329 mm	
	6.5 mm	9.0 mm	330 mm	334 mm	
	7.0 mm	9.6 mm	349 mm	352 mm	
	7.5 mm	10.2 mm	353 mm	356 mm	
104202, Murphy, 1 eye sizes: I.D. 3.5–10.0 mm	8.0 mm	10.8 mm	360 mm	364 mm	
	8.5 mm	11.3 mm	384 mm	388 mm	
	9.0 mm	11.9 mm	393 mm	398 mm	
104204, Murphy with FlexiSlip sizes: I.D. 5.0–9.0 mm	9.5 mm	12.4 mm	394 mm	399 mm	
	10.0 mm	12.9 mm	395 mm	399 mm	

\* Length without connector

# ARMoured TRACHEAL TUBES

Ref. 104004



**ARMoured TRACHEAL TUBE, STRAIGHT, MADE OF PVC, WITH LOW-PRESSURE CUFF**  
nasal/oral  
sizes: I.D. 5.0–9.0 mm

valve for Luer and Luer-lock syringes, radiopaque, blue pilot balloon, graduated

- latex-free
- sterile



## ARMoured TRACHEAL TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
104004 fixed white connector	5.0 mm	7.8 mm	17.5 mm	253.0 mm	2
	5.5 mm	8.3 mm	22.0 mm	253.5 mm	
	6.0 mm	8.8 mm	22.0 mm	313.5 mm	
	6.5 mm	9.3 mm	26.5 mm	313.5 mm	
	7.0 mm	9.8 mm	26.5 mm	314.0 mm	
	7.5 mm	10.3 mm	28.5 mm	314.0 mm	
	8.0 mm	10.8 mm	28.5 mm	314.0 mm	
	8.5 mm	11.3 mm	29.5 mm	354.0 mm	
	9.0 mm	11.8 mm	29.5 mm	354.5 mm	



104000  
blue mounted  
connector  
**FO intubation**

**Tip:**  
Ref. 104000 without connector can be used for fiberoptic intubation along with a Mainz universal adaptor (see page 48).

**Attention:**  
We recommend spraying the inside of the tracheal tube and the outside of the suction catheter with SILKOSPRAY (Ref. 556000, see page 49) prior to use.

**PAEDIATRIC ARMoured TRACHEAL TUBE, STRAIGHT, MADE OF PVC, WITHOUT CUFF**  
nasal/oral  
sizes: I.D. 2.5–6.5 mm

mounted connector, radiopaque, graduated

- latex-free
- sterile



## PAEDIATRIC ARMoured TRACHEAL TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	QTY
103600	2.5 mm	4.1 mm	184.0 mm	2
	3.0 mm	4.8 mm	184.5 mm	
	3.5 mm	5.5 mm	184.5 mm	
	4.0 mm	6.0 mm	194.5 mm	
	4.5 mm	6.5 mm	198.0 mm	
	5.0 mm	7.2 mm	208.0 mm	
	5.5 mm	7.7 mm	208.0 mm	
	6.0 mm	8.2 mm	208.0 mm	
	6.5 mm	8.9 mm	268.5 mm	



\* Length without connector

## SILKOCLEAR FLEX

**ARMOURED TRACHEAL TUBE  
MADE OF SILICONE, WITH HIGHLY  
FLEXIBLE THIN-WALLED CUFF**  
nasal/oral, 1 eye (Murphy)  
sizes: I.D. 3.5–10.0 mm

fixed connector, valve for Luer and  
Luer-lock syringes, black position  
indicator for correct tube placement,  
pilot balloon, smooth inner tube  
coating, graduated

- latex-free
- sterile



### SILKOCLEAR FLEX

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
105702	3.5 mm	5.0 mm	7.0 mm	226.0 mm	2
	4.0 mm	5.7 mm	7.0 mm	246.0 mm	
	4.5 mm	6.3 mm	7.0 mm	250.5 mm	
	5.0 mm	6.7 mm	8.5 mm	268.0 mm	
	5.5 mm	7.3 mm	12.0 mm	299.0 mm	
	6.0 mm	8.0 mm	12.0 mm	330.5 mm	
	6.5 mm	8.7 mm	12.0 mm	331.0 mm	
	7.0 mm	9.3 mm	13.5 mm	332.0 mm	
	7.5 mm	10.0 mm	13.5 mm	342.5 mm	
	8.0 mm	10.7 mm	16.0 mm	345.5 mm	
	8.5 mm	11.3 mm	16.0 mm	376.0 mm	
	9.0 mm	12.0 mm	16.0 mm	375.0 mm	
	9.5 mm	12.7 mm	18.0 mm	375.5 mm	
10.0 mm	13.3 mm	18.0 mm	378.0 mm		

**PAEDIATRIC ARMOURD  
TRACHEAL TUBE, STRAIGHT,  
MADE OF SILICONE,  
WITHOUT CUFF**  
nasal/oral  
sizes: I.D. 2.0–6.5 mm

fixed connector, cupped atraumatic  
tip, radiopaque, graduated

- latex-free
- sterile



### PAEDIATRIC ARMOURD TRACHEAL TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	QTY
105400	2.0 mm	3.3 mm	223.0 mm	2
	2.5 mm	3.9 mm	224.0 mm	
	3.0 mm	4.3 mm	225.5 mm	
	3.5 mm	5.0 mm	226.0 mm	
	4.0 mm	5.7 mm	246.0 mm	
	4.5 mm	6.3 mm	250.5 mm	
	5.0 mm	6.7 mm	263.0 mm	
	5.5 mm	7.3 mm	294.5 mm	
	6.0 mm	8.0 mm	326.0 mm	
	6.5 mm	8.7 mm	327.5 mm	

\* Length without connector

**Please note:** SILKOSPRAY is not suitable for devices made of silicone.  
If you require a lubricant, we recommend commercially available gels.

**ARMoured TRACHEAL TUBE,  
STRAIGHT, MADE OF SILKOLATEX,  
WITH SILICONE CUFF**

nasal/oral

sizes: I.D. 3.0–10.5 mm

O.D. in CH 18.0–42.0

mounted connector, radiopaque,  
graduated, with cuff, pilot balloon  
and universal injection port

- sterile

**Tip:** Ref. 103040 without connector can be used  
for fiberoptic intubation along with a Mainz  
universal adaptor (see page 58)

**Attention:**

These products contain natural rubber latex which  
may cause allergic reactions.



**ARMoured TRACHEAL TUBE**

**RÜSCH**

REF.	ORDER SIZE	I.D.	O.D.	LENGTH*	QTY
103040	Ch. 18	3.0 mm	6.0 mm	229.0 mm	5
	Ch. 20	3.5 mm	6.7 mm	232.5 mm	
	Ch. 22	4.0 mm	7.3 mm	232.5 mm	
	Ch. 24	5.0 mm	8.0 mm	253.0 mm	
	Ch. 26	5.5 mm	8.7 mm	253.0 mm	
	Ch. 28	6.0 mm	9.3 mm	313.5 mm	
	Ch. 30	6.5 mm	10.0 mm	313.5 mm	
	Ch. 32	7.0 mm	10.7 mm	314.0 mm	
	Ch. 34	8.0 mm	11.3 mm	354.0 mm	
	Ch. 36	8.5 mm	12.0 mm	354.5 mm	
	Ch. 38	9.0 mm	12.7 mm	354.5 mm	
	Ch. 40	9.5 mm	13.3 mm	355.0 mm	
Ch. 42	10.5 mm	14.0 mm	355.0 mm		

**BRONCHOFLEX-SET**

**ARMoured TRACHEAL/  
INSERTION TUBE MADE OF PVC  
FOR USE IN BRONCHOSCOPY**

(with separate oxygen tubing)

sizes: I.D. 7.5 mm and 8.5 mm

mounted connector, radiopaque,  
approx. 33 cm long, double lumen,  
graduated, siliconised, with or  
without cuff

- latex-free
- sterile

**Attention:**

We recommend spraying the inside of the  
bronchial tube and the outside of the suction  
catheter with SILKOSPRAY (Ref. 556000,  
see p. 59) prior to use.



**BRONCHOFLEX-SET**

**RÜSCH**

REF.	ORDER SIZE/I.D.	TOTAL O.D.	LENGTH*	QTY
104100	7.5 mm	13.8 mm	319 mm	1
<b>SET COMPONENTS</b> <ul style="list-style-type: none"> <li>• armoured tracheal/insertion tube made of PVC</li> <li>• oxygen tubing made of PVC, approx. 110 cm long</li> <li>• plastic bite block with fixation flange</li> <li>• neck band</li> </ul>				



**BRONCHOFLEX-SET**

**RÜSCH**

REF.	ORDER SIZE/I.D.	TOTAL O.D.	LENGTH*	QTY
104103	8.5 mm	11.7 mm	359 mm	1
<b>SET COMPONENTS</b> <ul style="list-style-type: none"> <li>• armoured tracheal/insertion tube made of PVC with high-volume low-pressure cuff</li> <li>• oxygen tubing made of PVC, approx. 110 cm long</li> <li>• valve for Luer-lock syringes</li> <li>• neckband</li> </ul>				

\* Length without connector

## LASER TUBE



**PREFORMED TRACHEAL TUBE  
WITH DEFINED LASER-RESISTANCE  
MADE OF WHITE SOFT RUBBER  
FOR LASER INTERVENTION IN  
THE LARYNGEAL SPACE**

oral/nasal (Magill), 1 eye (Murphy)  
sizes: I.D. 4.0–8.0 mm

fixed white connector, two-way stop-cocks with Luer and Luer-lock adapter and locking mechanism, cupped atraumatic tip, two cuffs inside one another, two pilot balloons, non-kinking shaft, Laserguard film 17 cm long consisting of Merocel® foam and corrugated silver foil

- sterile

**Attention:**

*These products contain natural rubber latex which may cause allergic reactions.*

### LASER TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø OUTSIDE	CUFF Ø INSIDE	LENGTH*	QTY
102004*	4.0 mm	8.0 mm	20 mm	20 mm	400 mm	2
	5.0 mm	9.0 mm	20 mm	20 mm	400 mm	
	6.0 mm	10.8 mm	24 mm	24 mm	400 mm	
	7.0 mm	12.3 mm	28 mm	28 mm	400 mm	
	8.0 mm	13.4 mm	28 mm	28 mm	400 mm	

\* Length without connector



## CLEAR MICROLARYNGEAL TUBE



### MADE OF PVC WITH HIGH-VOLUME LOW-PRESSURE CUFF

oral/nasal

sizes: I.D. 4.0–6.0 mm

thinner tube body designed for interventions in the laryngeal space

semi-seated connector, valve for Luer and Luer-lock syringes, continuous X-ray marker, blue pilot balloon, graduated

- latex-free
- sterile



### MICROLARYNGEAL TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
112460	4.0 mm	6.0 mm	31 mm	360 mm	10
	5.0 mm	7.3 mm	31 mm	360 mm	
	6.0 mm	8.7 mm	31 mm	360 mm	

\* Length without connector

## OXFORD TUBE

### CLEAR ANATOMICALLY SHAPED NON-KINKING TUBE MADE OF PVC WITH LOW-PRESSURE CUFF

sizes: I.D. 7.0–10.0 mm

semi-seated connector, valve for Luer and Luer-lock syringes, atraumatic tip, continuous X-ray marker, blue pilot balloon, graduated

- latex-free
- sterile



### OXFORD TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
112880	7.0 mm	10.0 mm	26 mm	229.5 mm	10
	8.0 mm	11.3 mm	28 mm	237.0 mm	
	9.0 mm	12.7 mm	29 mm	243.5 mm	
	10.0 mm	14.0 mm	29 mm	250.5 mm	



## SOFT RUBBER TRACHEAL TUBE



### MADE OF SOFT RED RUBBER, WITH CUFF

Magill  
sizes: I.D. 2.5–11.0 mm

semi-coated with SILKOLATEX®,  
cuff and pilot balloon made of  
SILKOLATEX®, universal injection  
port, graduated

- suitable for MRT
- without metal

**Attention:**

*These products contain natural rubber latex  
which may cause allergic reactions*



### SOFT RUBBER TRACHEAL TUBE

**RÜSCH**

REF.	ORDER SIZE/I.D.	O.D.	LENGTH*	QTY
102000	2.5 mm	4.0 mm	140 mm	10
	3.0 mm	4.7 mm	160 mm	
	3.5 mm	5.3 mm	180 mm	
	4.0 mm	6.0 mm	200 mm	
	4.5 mm	6.7 mm	220 mm	
	5.0 mm	7.3 mm	240 mm	
	5.5 mm	8.0 mm	270 mm	
	6.0 mm	8.7 mm	280 mm	
	6.5 mm	9.3 mm	290 mm	
	7.0 mm	10.0 mm	300 mm	
	7.5 mm	10.7 mm	310 mm	
	8.0 mm	11.3 mm	320 mm	
	8.5 mm	12.0 mm	330 mm	
	9.0 mm	12.7 mm	340 mm	
	9.5 mm	13.3 mm	350 mm	
	10.0 mm	14.0 mm	360 mm	
10.5 mm	14.7 mm	360 mm		
11.0 mm	15.3 mm	360 mm		

\* Length without connector



### STANDARD CONNECTOR

**RÜSCH**

REF.	ORDER SIZE/I.D.	DESCRIPTION	QTY
501003	1.5 mm–11.0 mm	made of plastic, latex-free, 15 mm connection (O.D.)	10



# AID INSTILLATION CATHETER

**FOR DRUG INSTILLATION DURING CPR AND SURFACTANT APPLICATION IN NEONATOLOGY**

with sealable Luer-lock adapter, multiple perforations over a length of 1 cm, opening 1 cm from distal end of the catheter, with intubation stylet, graduated

- latex-free
- sterile

# EDGAR TUBE

**CLEAR TRACHEAL TUBE MADE OF PVC FOR ENDO-BRONCHIAL DRUG & GAS APPLICATION DURING RESUSCITATION**  
nasal/oral (Magill)

semi-seated connector, additional instillation channel with sealable Luer-lock adapter, continuous X-ray marker, graduated

**Ref. 111480**

sizes: I.D. 6.5–10.0 mm with low-pressure cuff, valve for Luer and Luer-lock syringes and black position indicator for correct tube placement


**Ref. 111380**

sizes: I.D. 2.5–6.0 mm without cuff, with black tip

- latex-free
- sterile



AID ADAPTER SET			RÜSCH
REF.	SIZE/O.D.	DESCRIPTION	QTY
111100	Ch. 4	catheter with angled connector	1
		<b>SET COMPONENTS</b> • AID instillation catheter • angled plastic connector	

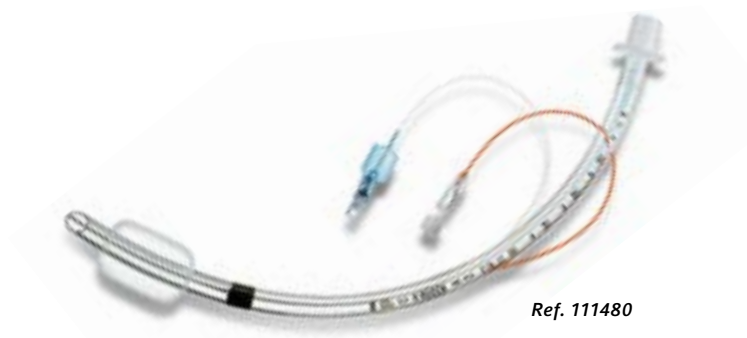

EDGAR TUBE					RÜSCH
REF.	ORDER SIZE/I.D.	O.D.	CUFF Ø	LENGTH*	QTY
111480	6.5 mm	9.3 mm	24 mm	290 mm	10
	7.0 mm	10.0 mm	26 mm	300 mm	
	7.5 mm	10.7 mm	26 mm	310 mm	
	8.0 mm	11.3 mm	28 mm	320 mm	
	8.5 mm	12.0 mm	28 mm	330 mm	
	9.0 mm	12.7 mm	29 mm	340 mm	
	9.5 mm	13.3 mm	29 mm	350 mm	
	10.0 mm	14.0 mm	29 mm	360 mm	
111380	2.5 mm	4.0 mm	–	140 mm	10
	3.0 mm	4.7 mm	–	160 mm	
	3.5 mm	5.3 mm	–	180 mm	
	4.0 mm	6.0 mm	–	200 mm	
	4.5 mm	6.7 mm	–	220 mm	
	5.0 mm	7.3 mm	–	240 mm	
	5.5 mm	8.0 mm	–	270 mm	
	6.0 mm	8.7 mm	–	280 mm	

\* Length without connector